


Thank you for purchasing this product. This manual describes the operating instructions and safety precautions for using this device. Please be sure to read this manual carefully before using the device. After reading, keep the manual in a safe place for future reference.

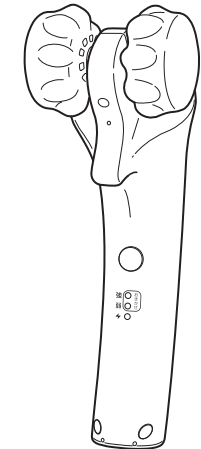

1st edition, September 2017

Safety Precautions (Be sure to always follow these precautions)

In order to prevent personal injury or damage to property involving you or any other people, please be sure to observe the precautions outlined below.

■ The following symbols show the degree of injury or damage which may result through incorrect use or failure to observe the precautions in the text.

| | | |
|--|----------------|---|
| | Warning | This symbol indicates the possibility of death or serious injury. |
| | Caution | This symbol indicates the possibility of bodily harm or damage to property. |
| | | This symbol indicates a "Prohibited" action which must not be performed. |
| | | This symbol indicates a "Mandatory" action which must be performed. |

| | | |
|--|----------------|---|
| | Warning | <p>Any person whose condition corresponds to any of the following should consult a doctor prior to use.</p> <ul style="list-style-type: none">• People using an implanted medical electric device such as a pacemaker that is susceptible to electromagnetic interference• People receiving treatment from a physician, or people who feel a particular abnormality in their body• People with malignant tumors• People with heart problems• Pregnant women or those who have just given birth• People with impairments in perception arising from advanced peripheral circulatory disorders due to diabetes.• People with wounds on their skin• People who require rest in bed• People with body temperatures of 38 ° C or higher (at fever stage)* Example 1: Acute inflammatory symptoms (such as fatigue, chills, blood pressure fluctuations)* Example 2: Debilitated condition• People with osteoporosis or those with acute (painful) disorders, such as spinal fractures, sprains, pulled muscles• People with an abnormal backbone or those whose backbone is bent to the right or left• People with a herniated disc• People feeling a physical abnormality or those under treatment at medical institutions• People who feel no effects despite using the device for a while <p>Such use could cause an accident or poor physical health.</p> |
|--|----------------|---|

2

| | | |
|--|----------------|---|
| | Warning | <p>Any person whose condition corresponds to any of the following should consult a doctor prior to use.</p> <ul style="list-style-type: none">• People whose doctors have prohibited them from massages* Example: People with thrombosis (embolism), serious aneurysm, acute varicose veins, various types of dermatitis or skin infection (including inflammation of subcutaneous tissues) <p>Such use could worsen symptoms.</p> |
| | | <p>Do not use this device for any purpose other than the intended one. Do not use this device on your head, face, breast, belly, bony parts (such as knee or elbow) or genitals.</p> <p>Such use could cause an accident or injury.</p> |
| | | <p>Do not press the roller or pressing point too hard during use.</p> <p>Such use could cause an accident or injury.</p> |
| | | <p>People who are unable to express themselves or use the device by themselves should not be allowed to use the device by themselves. Do not allow children to use or play with this device. Keep this device out of reach of small children when using or storing the device.</p> <p>Otherwise, an accident, injury or device failure may be caused. If a person has swallowed a part, consult a physician immediately.</p> |
| | | <p>Do not use the device by standing or sitting on the main unit.</p> <p>Otherwise, an accident, injury or device failure may be caused.</p> |
| | | <p>Do not use this device simultaneously with any other therapeutic apparatus.</p> <p>Such use could cause discomfort or poor health.</p> |
| | | <p>Do not disassemble, repair or modify this device.</p> <p>Fire or electric shock may be caused.</p> |

Caution

| | |
|--|--|
| | <p>Do not use this device for longer than ten minutes in a day. Do not continuously use the device for longer than three minutes in any area.</p> <p>Such use could cause an accident or poor physical health.</p> |
| | <p>Always hold the device with your hand when using the device.</p> <p>Otherwise, an accident, injury or device failure may be caused.</p> |

3

| | | |
|--|----------------|---|
| | Caution | |
| | | <p>If you experience any unusual physical sensations in your body during use, stop immediately and consult a physician.</p> <p>Such use could cause an accident or poor health.</p> |
| | | <p>After use, stop operation with the power button.</p> <p>Such use could cause an accident, injury, damage of the product due to the falling of the product.</p> |
| | | <p>If a device failure occurs, stop immediately, turn off the power, and contact the dealer.</p> <p>Such use could cause an accident or injury.</p> |
| | | <p>Do not use this device when:</p> <ul style="list-style-type: none">• You are very tired;• You are drunk;• You have eczema. <p>Such use could cause poor health.</p> |
| | | <p>Do not use this device within one hour before or after a meal.</p> <p>Such use could cause poor health.</p> |
| | | <p>Do not sleep during use.</p> <p>Such use could cause an accident or injury.</p> |
| | | <p>Do not use this device in bed or under a blanket. Do not use this device near a fire such as a stove or while smoking.</p> <p>Such use could cause burning, overheating of the motor or a fire hazard.</p> |
| | | <p>Do not use this device when you are wearing hard objects such as jewelry like necklaces or watches.</p> <p>Such use may cause injury.</p> |
| | | <p>Do not drop this device nor subject it to strong impacts.</p> <p>Such use may cause an electric shock, injury or device failure.</p> |
| | | <p>Do not use thinner, benzene, alcohol or hot water to clean the device.</p> <p>Discoloration or deformation may be caused.</p> |
| | | <p>Do not use the device in a wet or humid place such as the bathroom or in water.</p> <p>Such use may cause electric shock, device failure or fire.</p> |
| | | <p>Do not use the device while driving a car, riding a bicycle or motorcycle.</p> <p>Such use could cause an accident or problem.</p> |

4

| | | |
|--|----------------|---|
| | Caution | |
| | | <p>Avoid touching the DC connector socket with any part of your body such as a hand or a finger.</p> <p>Otherwise, electric shock or device failure may be caused.</p> |
| | | <p>Do not insert any metallic object in the DC connector socket.</p> <p>Otherwise, overheating or fire may be caused.</p> |
| | | <p>Avoid excessive pulling or bending of the dedicated USB cable, or placing heavy objects on the cable.</p> <p>Otherwise, cable damage or breakage, fire or electric shock may be caused.</p> |
| | | <p>Do not use the device if the dedicated USB cable is damaged, or the DC connector socket is loose.</p> <p>Otherwise, fire or electric shock may be caused.</p> |
| | | <p>Do not use a travel converter or any other power converter to charge the device.</p> <p>Otherwise, overheating or fire may be caused.</p> |
| | | <p>Do not charge in a car using the cigarette lighter receptacle or any other means.</p> <p>Otherwise, fire, overheating or device failure may be caused.</p> |
| | | <p>Do not touch the cable with a wet hand.</p> <p>Otherwise, electric shock, device failure or accident may be caused.</p> |
| | | <p>Do not use the supplied dedicated USB cable for any purpose other than charging the main unit.</p> |
| | | <p>Do not charge this device when the main unit is wet.</p> <p>Such use may cause overheating, fire or electric shock.</p> |
| | | <p>Please confirm that the power button works properly before use.</p> |
| | | <p>When using for the first time, or when the device has been out of use for a long period of time, follow the instruction manual and confirm it operates normally and safely prior to use.</p> |
| | | <p>Before use, make sure that there is no crack in the main unit and that the power button is not broken.</p> <p>Cracks and breakage affect waterproofing performance and may cause electric shock, device failure or fire.</p> |

5


| | | |
|--|----------------|--|
| | Caution | |
| | | <p>Make sure the dedicated USB cable is connected securely when charging.</p> <p>If connected incompletely while charging, fire or an electric shock may be caused.</p> |
| | | <p>Be sure to hold the connector to pull out the dedicated USB cable from the main unit.</p> <p>Do not pull the cable. Otherwise, electric shock, short circuit or device failure may be caused.</p> |
| | | <p>Be sure to charge the battery at least once every three months.</p> <p>If the device is fully discharged, it may not be rechargeable.</p> |
| | | <p>If a power failure has occurred while charging, immediately stop charging.</p> <p>Otherwise, electric shock, device failure or accident may be caused.</p> |
| | | <p>When disposing of this product, please do so in accordance with the regulations of your municipality.</p> |
| | | <p>Carefully read this instruction manual before use. Make sure there is no problem with this device before use.</p> |
| | | <p>After use, clean the device for the next use and store it in a low humidity place.</p> |

6

Package contents

Please make sure all components are included. We aim to ensure that all our products are carefully packaged. Should any item(s) be missing or damaged, immediately contact the dealer you purchased the product from.

Main Unit

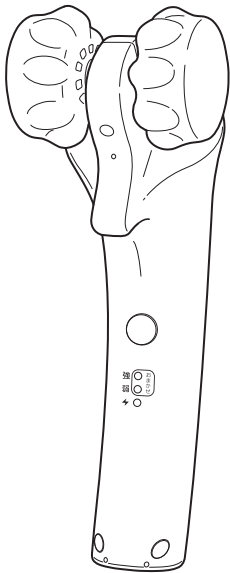

Dedicated USB cable




Instruction manual (with warranty) - this document

7

Part names

Main Unit

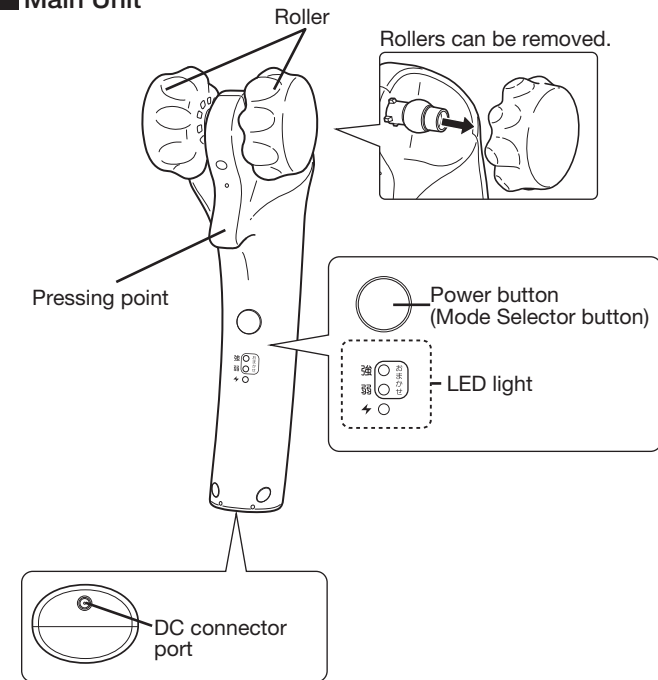

8

Dedicated USB cable



8

Preparation


Be sure to charge the battery fully on first use (following purchase or after a long period without use) or if the batteries have run out during use.

Charging


Charging time: approximately 2 hours

- Connect the DC connector of the dedicated USB cable to the DC socket on the main unit.
 - Connect Type-A USB connector of the dedicated USB cable to any of the USB ports shown below to charge the battery.
- * When correctly connected, the LED indicator on the main unit will blink as shown on the lower right.

Main Unit

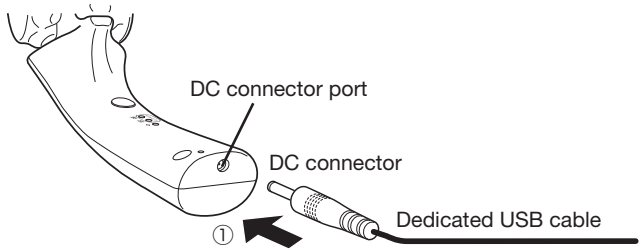

- Connecting to a PC
- USB 2.0/3.0 port of a PC conforming to the IEC60950 standard

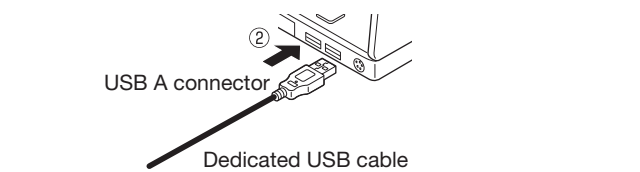

9

- When fully charged, the LED indicator will stop blinking and stay lit.
- When charging has finished, disconnect the dedicated USB cable from the main unit, following the steps backwards from ② to ① for safety. Once the battery is fully charged, it can be used for approximately 3 hours.
- * This will change depending on the type of operation and usage environment.

CAUTION :

- Use of travel converters or other power converters for charging may cause overheating or fire.
- Do not charge the device when the DC connector socket of the main unit is wet.
- The device does not operate while charging.
- The operating time of a fully charged battery becomes shorter with repeated charging of the main unit. If the operating time becomes extremely short, the battery has reached the end of its life. However, the battery is not replaceable.

LED light


The LED indicator status changes as shown below depending on the remaining battery power.

| | |
|--|---|
| [0 to 30%] 強 (off) 弱 (off) | [30 to 80%] 強 (off) 弱 (off) |
| [80% or higher] 強 (off) 弱 (off) | [Charging complete] 強 (off) 弱 (off) |
| <div>○ Light is off</div> <div>● Light is on</div> | |


10

How to use

Press and roll the rollers over the entire area where you wish to apply the device.

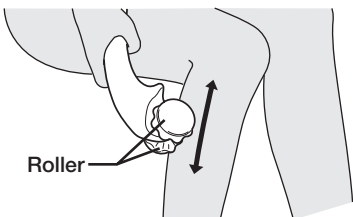


Press the pressing point onto the area where you wish to apply the device.

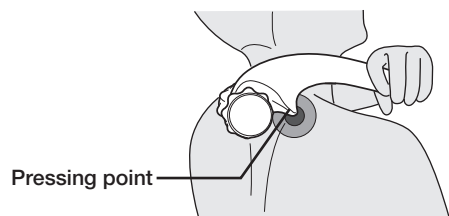


Turn on the vibrator, then roll the rollers or press the pressing point.

- Turn on the power.
- Press and hold the power button for about one second. The device will turn on and the vibrator will start. The Weak mode LED light turns on.

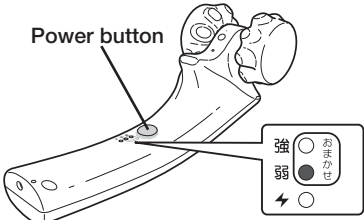


11

- Select mode.
- Press the power button to select a mode during operation. When the power button is pressed, the modes change in the sequence shown below.
- * The device always starts in the Weak mode when it is turned on.

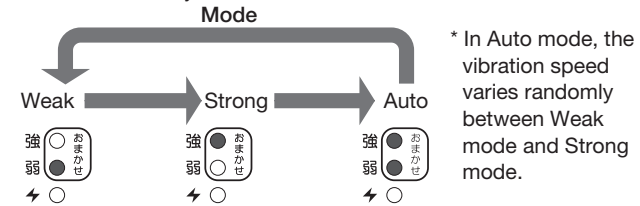

* In Auto mode, the vibration speed varies randomly between Weak mode and Strong mode.

- Turn off the power.
- The device turns off automatically about five minutes after the start of operation.
- * Press and hold the power button for about one second to turn off the device during operation.

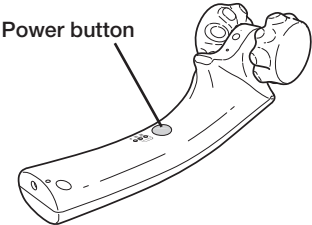

CAUTION :

- If you feel pain during use, apply a towel or the like underneath the device.
- Do not leave the device in high temperature or humidity areas, under direct sunlight, nor within reach of young children.

Recommended time of use

When your body is warm, such as after a bath. Since this product has waterproofing performance equivalent to JIS IPX6, it can be safely used on a wet body.

12

