Manual(English) ▶▶▶ p.2 - p.24

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
For residential use

Manual (ENGLISH)

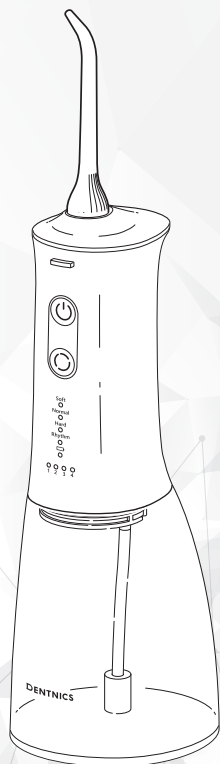
Read this manual and use the product safely and correctly.
After reading, keep the manual in an easily accessible location for future reference.

*In this manual, HC-JF03WH/HC-JF03BK will be referred to as “the product” except in some places.

Package contents

The product package includes the items below. Check if they are all supplied. Although every effort has been made to ensure proper packaging, contact the retailer where this product was purchased if any shortage or breakage is discovered.

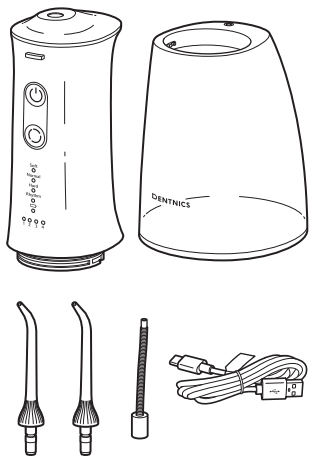- ☐ Main unit × 1
- ☐ Flosser nozzle × 2
- ☐ Water tank × 1
- ☐ Dedicated USB charging cable × 1
- ☐ Water hose × 1

*The water hose is pre-connected to the main unit.

Table of contents

Package contents	1	Cleaning and storage.....	13
Safety precautions	2	After each use.....	13
Names and functions of components .	5	Cleaning the product (weekly or so).....	14
Before use	7	Cleaning the water hose.....	15
Charging the product	7	Storing the product	15
How to use.....	9	Lithium-ion battery	16
Setting the flosser nozzle	9	Troubleshooting.....	17
Pouring water.....	9	FAQ	19
Selecting a mode	10	Replacement parts.....	20
Setting the water jet level	10	Specifications.....	20
Flossing your mouth.....	11	Regular inspections	21
		Compliance status	22

Safety precautions

Read the safety precautions carefully and use the product correctly. The precautions in this section describe the actions and behaviors to take to prevent hurting yourself and others or damaging property. The meaning of the symbols used in this manual is explained below.

● Meaning of symbols

 **WARNING** indicates matters that have a “probability of death or serious injury can be anticipated”.

 **CAUTION** indicates matters that have a “probability of injury or property damage can be anticipated”.

 General precautions

 Do not touch with wet hands

 Do not wet

 General prohibitions

 Do not disassemble

 Follow instructions

 indicates “matters that require careful attention”.

 indicates “recommendations” and “reference information”.

WARNING

 **Do not use the product for any purpose not described in the manual.**

Doing so may cause a fire, ignition, electric shock, burn, or injury.

 **Do not disassemble, repair, or modify the product.**

Doing so may cause a malfunction, fire, electric shock, or injury. In the event of a malfunction, stop using the product and disconnect the dedicated USB charging cable. Next, consult the retailer where this product was purchased.

 **Do not put the product in a fire, heat it, or leave it in an extremely hot place.**

The product contains a rechargeable battery. Being extremely hot, it may cause heat generation, ignition, rupture, or injury.

 **Do not use the product near medical or home therapy equipment.**

Do not use the product near implanted medical devices (e.g., a cardiac pacemaker and a defibrillator), life support and management devices (e.g., a heart-lung machine and a respirator), or wearable medical devices (e.g., a wearable automatic defibrillator and an ambulatory infusion pump). Doing so may cause the devices to malfunction.

 **Do not allow children or those who are incapable of, requiring assistance in, or unfamiliar with operating the product to use it alone.**

Doing so may cause an electric shock, burn, or injury.

 **Do not allow water to spill over the dedicated USB charging cable or USB AC adapter (commercially available). Do not fully wash them.**

Doing so may cause a malfunction, fire, or electric shock.

 **Do not use or store the dedicated USB charging cable and USB AC adapter (commercially available) in a bathroom or any other humid place.**

Doing so may cause a fault current or short circuit, resulting in a fire or electric shock.

Safety precautions (continued)

WARNING

-
-  **Always use the dedicated USB charging cable.**
Failure to do so may cause a malfunction, error, fire, ignition, electric shock, or burn.
-
-  **Do not allow conductive foreign objects (such as metal pieces, pencil lead, or foreign particles) to get inside the USB terminal.**
They may cause a malfunction, fire, ignition, electric shock, or burn.
-
-  **Immediately stop using the product in the event of an abnormality or malfunction.**
Failure to do so may cause a fire, ignition, electric shock, burn, or injury.
- There is a strange odor or noise.
 - The dedicated USB charging cable or USB AC adapter (commercially available) gets abnormally hot.
 - You feel tingling when touching the main unit.
 - The main unit has cracked or broken.
- If you find symptoms above, immediately stop using the product and disconnect the dedicated USB charging cable. Next, consult the retailer where this product was purchased.
-
-  **Do not use wet hands when connecting or disconnecting the dedicated USB charging cable and USB AC adapter (commercially available).**
Doing so may cause a fire, electric shock, or burn.
-
-  **Always use a USB AC adapter (commercially available) that is USB Type-A with a rated output of 5 V/1 A or higher.**
Failure to do so may cause a malfunction, error, fire, ignition, electric shock, or burn.
-
-  **Do not use a dedicated USB charging cable that is damaged or has a loose USB terminal.**
Doing so may cause a fire, ignition, electric shock, or burn.
-
-  **Do not use the dedicated USB charging cable in such a manner that it may generate heat.**
Doing so may cause a fire, ignition, electric shock, or burn.
- Use it while coiled up
 - Modify or damage it
 - Apply excessive force to it by pulling, bending strongly, or twisting it
 - Place a heavy object on it
 - Leave it squeezed between objects or close to a high-temperature object
-
-  **Do not wind the dedicated USB charging cable around the main unit.**
Doing so may cause a broken cable, fire, ignition, electric shock, or burn.
-
-  **Do not use the product or put accessories or parts within reach of young children.**
Doing so may cause an electric shock, injury, or malfunction.
-

Safety precautions (continued)

CAUTION

-  **Do not use the product on unstable or slippery surfaces.**
Doing so may cause an injury or malfunction.
-  **Do not drop, step on, throw, or subject this product to strong forces, shocks, or vibrations.**
Doing so may cause a broken dedicated USB charging cable or USB terminal, resulting in a fire, ignition, electric shock, or burn.
-  **Disconnect the dedicated USB charging cable and unplug the USB AC adapter (commercially available) from the outlet except while charging the product.**
Failure to do so may cause a malfunction, fire, or electric shock.
-  **Use the product for residential use.**
The product is designed for residential use.
-  **Do not use the product for any purpose not described in the manual, or for young children or pets.**
Doing so may cause a malfunction, fire, or burn.
-  **Do not use the product for any purpose other than flossing your mouth.**
Water jets may cause an injury when applied to your eyes, nose, ears, throat, under the tongue, or any other delicate part.
-  **If you are suspected of gum disease or are concerned about the health of your teeth and mouth, consult a dentist before use.**
Failure to do so may cause injury or damage to your teeth and gums.
-  **Do not use the product if you feel numb in your mouth.**
Doing so may cause injury or damage to your teeth and gums.
-  **Do not allow any part other than the tip of the flosser nozzle to come into contact with your teeth and gums.**
Doing so may cause injury or damage to your teeth and gums.
-  **Do not use hot water above 40 °C.**
Doing so may cause a burn.
-  **Do not use any liquid other than cold or lukewarm water, such as mouthwash and dental rinse.**
Doing so may cause a malfunction.
-  **Do not share a flosser nozzle with other people including your family members.**
Doing so may cause an infection or inflammation through saliva and mucosa.
-  **After use, clean and dry the main unit, water tank, flosser nozzle, and water hose.**
Failure to do so may allow mold and bacteria to grow, which can be harmful to your health.

To dispose of the product and package, follow the local rules in your area.

Names and functions of components

● Main unit

Mode selector button

Pressing the button shows the selected mode. Each time you press it, the mode switches in order of soft mode, normal mode, hard mode, and rhythm mode. You can confirm the selected mode with the light of the mode indicator. Since the last mode used is stored in the memory, it resumes the next time the power is turned on. If the last mode used is hard mode, the mode switches in order of hard mode, rhythm mode, soft mode, and normal mode.

Charging indicator (red)

Blinks when the battery level is low, blinks when charging is in progress, and stays lit when charging is completed.

Water jet level (blue)

Pressing and holding the mode selector button for two seconds during flossing switches the water jet level in order of 1, 2, 3, and 4. Level 1 is the lowest and Level 4 is the highest.

Flosser nozzle attachment part

Flosser nozzle release button

Power button (On/Off)

Mode indicator (blue)

Soft mode

Gentle water jets treat teeth and gums. Optimal when you use the product for the first time or have sensitive gums.

Normal mode

For those who are accustomed to soft mode. Moderate water jets wash off plaque on the surfaces and food particles between teeth.

Hard mode

Strong water jets wash off plaque on the surfaces and food particles between teeth.

Rhythm mode

Soft mode and hard mode alternate to stimulate and treat gums.

Lid of USB charging port

Water hose

Water tank

USB charging port

Connect the dedicated USB charging cable (USB Type-C® connector) supplied.

Lock/unlock

Turn the water tank to  to remove it from the main unit. Turn it to  to lock it.

Water port


Lid of water port

Names and functions of components (continued)

● Flosser nozzle x 2


*Flosser nozzles are consumables.

- Refer to “Replacement parts ⇒ p. 20”.

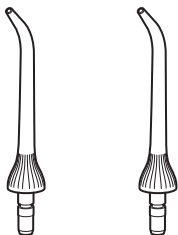● Dedicated USB charging cable

USB Type-C® connector

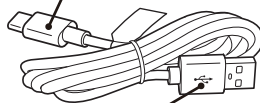USB Type-A connector

*This product does not include a USB AC adapter. Use a commercially available USB AC adapter that is USB Type-A with a rated output of 5 V/1 A or higher.

Before use

[Cautions for use]

- Use tap water or lukewarm water (40 °C or lower). Do not use salt water, mouthwash, dental rinse, or other solutions and chemicals for mouth cleansing.
- Do not use a dryer or fan heater to dry the main unit and accessories.
- Store the product in a less humid place after use. Do not leave it in a humid place such as a bathroom.

[Precautions before use]

● About this product

- Use the product in combination with tooth brushing. Plaque building up on the tooth surfaces cannot be removed by water jets alone, so tooth brushing is essential. It is recommended that you brush your teeth before flossing your mouth with the product.

*Effects may vary from person to person.

● Bleeding from gums

- With unusual stimulus on gums, you may experience bleeding for the first few weeks after starting use, even if you have healthy gums. If bleeding is intense or continues for longer than two weeks, it may be a symptom of abnormal gums (such as gum disease). If so, stop use immediately and consult a dentist.

● Precautions about waterproofness

- The main unit is waterproof (IPX7*) enough to allow use in a bathroom, but not completely waterproof. Do not use it under water.

*The product has passed the watertight test that evaluates its capability to prevent water entry that may cause adverse effects on the main unit by gradually immersing it one meter under water and leaving it for 30 minutes.

- The dedicated USB charging cable is not waterproof. Do not charge the product in a bathroom, near a humidifier, or in any other place that may be exposed to water.

■ Charging the product


- If you are using the product for the first time (just after purchase) or have left it unused for more than six months, be sure to charge it fully up to maximum.
- If the battery becomes low during use and the charging indicator blinks red, be sure to charge it fully up to maximum.
- This product does not include a USB AC adapter. Use a commercially available USB AC adapter that is USB Type-A with a rated output of 5 V/1 A or higher.

1

Open the lid of the USB charging port of the main unit and connect the dedicated USB charging cable to the USB charging port.

2

Connect the dedicated USB charging cable to a USB AC adapter (commercially available), and then to an outlet.

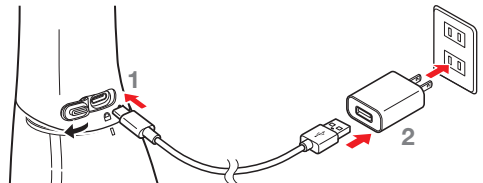
Wipe the USB charging port with a cloth and dry it well if it is wet. Leaving it wet or damp may cause a fault current.

Before use (continued)

3

The charging indicator starts blinking red slowly and charging starts.

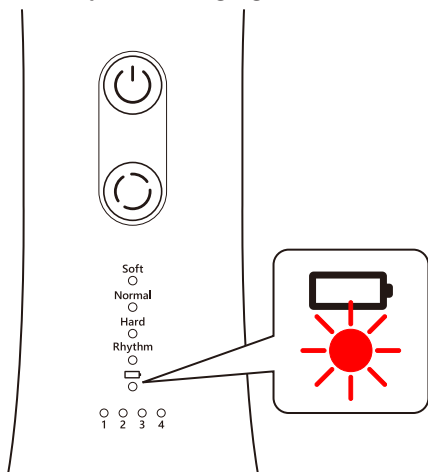4

If the charging indicator stays lit blue and charging is completed, disconnect the dedicated USB charging cable from the main unit.

[Charging indicator status]

Charging in progress:

The indicator blinks red slowly.

- It will take approximately three hours to fully charge a completely exhausted battery.

*The product cannot be used during charging.

Charging completed:

The indicator stays lit blue.

- With the battery fully charged, the maximum continuous operating time is approximately 40 minutes (in hard mode).
- *It depends on the ambient temperature and the use conditions.

Low battery: The indicator blinks red.

- The lamp blinking red when the power is turned on means it is time to charge.

- You can continue using the product for a while after the charging indicator starts blinking red, but the power will go out soon. Charging it early after blinking starts is recommended.
- If you are using the product for the first time (just after purchase) or have left it unused for more than six months, the charging time may change.
- Repeatedly charging the battery gradually shortens the operable time after a full charge. If the operable time has become extremely short, the battery has reached the end of its life. However, the battery cannot be replaced. Refer to “[Lithium-ion battery](#) ⇒ p. 16” for details.
- The product may not be charged or take more time to be charged outside the recommended operating temperature (5 °C to 40 °C).

- Immediately stop using the dedicated USB charging cable if it is broken. Purchase a USB Type-C — USB Type-A cable (commercially available) for replacement.

5

Be sure to close the lid of the USB charging port tightly to prevent water from entering the USB charging port.


How to use


Check the flosser nozzle for damage and deformation before use.

■ Setting the flosser nozzle

Hold the main unit and attach the flosser nozzle by pushing it until it clicks.

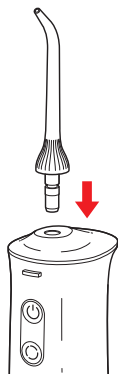
- Set the flosser nozzle securely so that its tip faces the front.

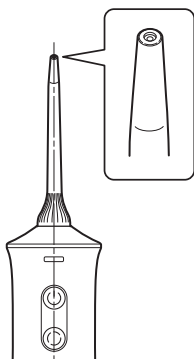■ Pouring water


1

Open the lid of the water port of the water tank.

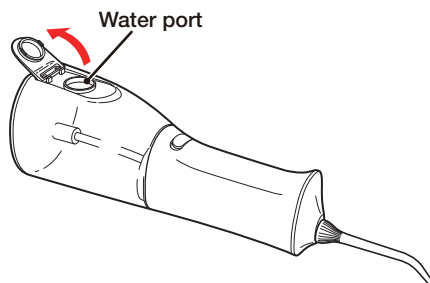2

Pour water into the water tank from the water port.

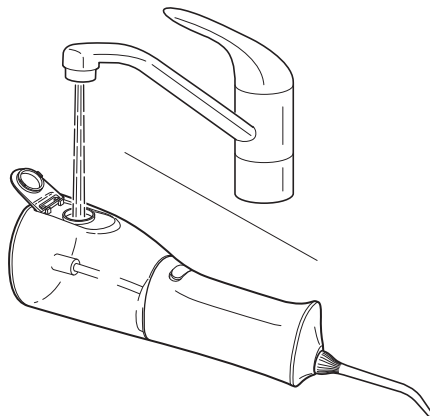- Always use fresh tap water or lukewarm water.
- If the water is too cold for you, use lukewarm water of 40 °C or lower.

3

Close the lid of the water port.

How to use (continued)

[Rinsing]

- If you use the product for the first time (just after purchase) or have left it unused for more than six months, rinse it in hard mode for approximately 10 seconds.

The recommended operating temperature is 5 °C to 40 °C. The product may not work outside the recommended operating temperature.

■ Selecting a mode

Press the mode selector button to select a mode.



- Pressing the mode selector button causes the mode indicator of the selected mode to light blue.
- To see the selected mode later, press the mode selector button again. The indicator of the selected mode lights blue. Then, it goes out in 10 seconds if left unoperated.
- Refer to [“Names and functions of components ⇒ p. 5”](#) for the details of the modes.
- “Soft mode” is recommended when you use the product for the first time.
- The selected mode is stored after the power is turned off. The last mode used resumes the next time the power is turned on. (Memory function)

■ Setting the water jet level

Press and hold the mode selector button for two seconds to set the water jet level.

- Pressing the mode selector button causes the water jet level indicator to light blue, indicating the set water jet level.
- To see the set water jet level later, press the mode selector button again. The water jet level indicator lights blue, indicating the set water jet level. Then, it goes out in 10 seconds if left unoperated.
- Starting with “Level 1” is recommended.
- The selected water jet level is stored after the power is turned off. The last water jet level used resumes the next time the power is turned on. (Memory function)

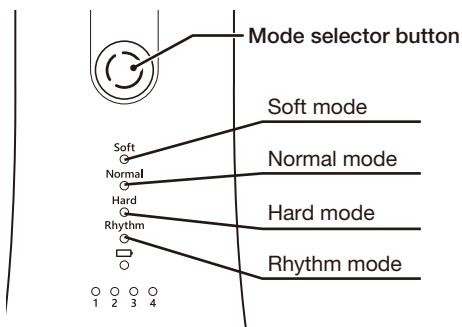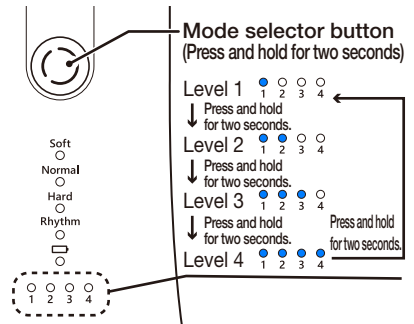How to use (continued)

■ Flossing your mouth

1

Above a basin, put the flosser nozzle into your mouth.


You can also use the product in a bathroom. Refer to "[Precautions before use] ⇒ p. 7" for the waterproofness.

2

Slightly close your mouth.

Water will splash if your mouth is open. Slightly close your mouth so that water can naturally flow down into the basin.

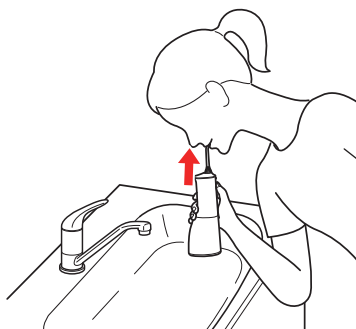3

Press the power button to start flossing.

- Press the power button once to start and once again to stop water jets.
- If the product is not operated for approximately two minutes after the power button is pressed once or the mode is switched, water jets automatically stop.
- You can also switch the mode by pressing the mode selector button during flossing.
- You can change the water jet level by pressing and holding the mode selector button for two seconds during flossing.
- Pour water into the water tank if it has become empty on the way.

How to use (continued)


[Tips for flossing]

- Change the mode and how to apply and move the flosser nozzle according to your purpose and the part to floss.

*If you are concerned about your gum condition, use soft mode.


To floss the gaps between teeth (recommended mode: hard mode)

- Keep the flosser nozzle over each gap for a short time during flossing.
- Apply water jets to the upper and lower teeth from their front and back.

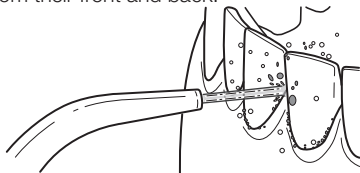
To floss the periodontal pockets (recommended mode: hard mode, normal mode, and soft mode)

- Apply water jets to the gumlines.
- Slowly move the flosser nozzle along the gums.
- Floss the front and back of the upper and lower teeth and the back of the molars.

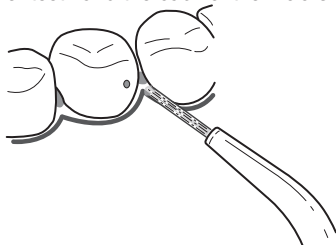To care for the gums (recommended mode: rhythm mode)

- Apply water jets to the gums to stimulate them.

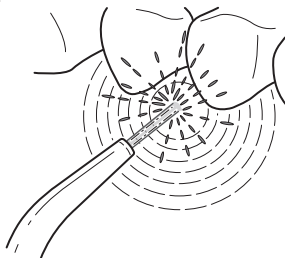4

After flossing, press the power button to turn the power off while the flosser nozzle is still in your mouth.

To prevent water splashing, wait until water feeding stops before taking the flosser nozzle out of your mouth.

If the product is not operated for approximately two minutes during continuous use, water jets automatically stop. After they automatically stop, press the power button again to resume. Keep the continuous use within approximately four minutes. After the continuous use, wait approximately 30 minutes before the next use.

Cleaning and storage

■ After each use

- After each use, be sure to throw away the water and clean the product to keep it hygienic.

1


After use, rinse the product in hard mode for approximately 10 seconds.

- Refer to “Selecting a mode ⇒ p. 10” for how to switch the mode.
- Refer to “[Rinsing] ⇒ p. 10” for details of rinsing.


2

Press the flosser nozzle release button to remove the flosser nozzle from the main unit.

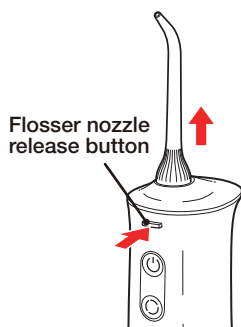3


Remove the water tank from the main unit.

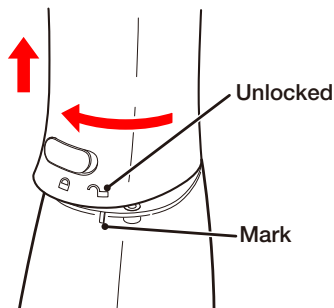- To remove the water tank, turn the main unit until the  (unlocked) mark on the main unit is aligned with the  mark on the water tank.
- Be careful not to touch the power switch of the main unit.

4

Empty the water tank.

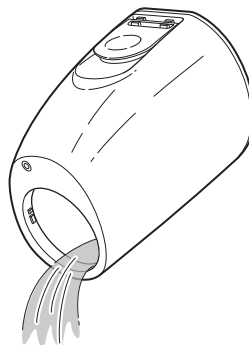Cleaning and storage (continued)

■ Cleaning the product (weekly or so)

- Do not use paint thinner, benzene, or alcohol. Doing so may cause a malfunction or may damage and discolor the parts.
- Do not press the power button during cleaning the product. Doing so may cause an electric shock, ignition, or malfunction.
- Do not use hot water of 50 °C or more on the main unit.
- Use neutral kitchen detergent if necessary, and rinse it off completely.


● Flosser nozzle, main unit, and water tank

1

Wash the flosser nozzle, main unit, and water tank with tap water.


If dirt is stubborn, wash them using neutral detergent diluted with water, and rinse it off completely.

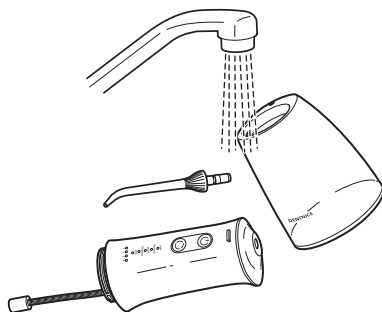Do not leave them under water.

2


Shake the main unit lightly to drain water from the inside.

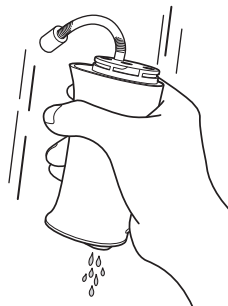Do not hit the main unit on the basin during draining.

3

Wipe water off with a soft, dry cloth and leave the parts to dry in a well-ventilated shade.

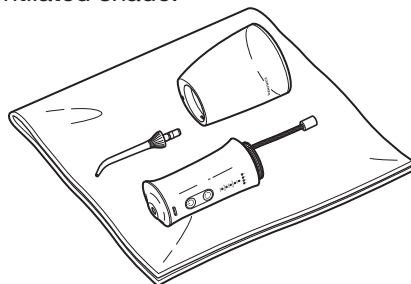Access cleaning details
from here →

<https://www.elecom.co.jp/r/n760>


Cleaning and storage (continued)

■ Cleaning the water hose

- You do not have to remove and clean the water hose every day.


If the water jet volume is low or foreign materials are adhering to the tip of the water hose, gently rub them off with a soft brush.

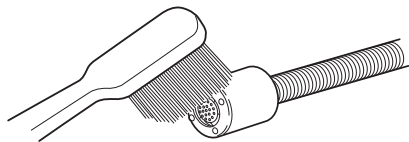
● To replace the water hose

- Replace the water hose in approximately six months after starting use for good hygiene.
- Refer to “Replacement parts ⇒ p. 20”.


To disconnect the water hose

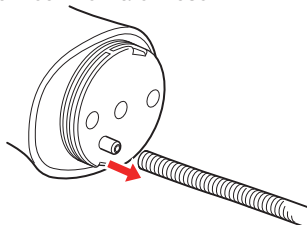To connect the water hose

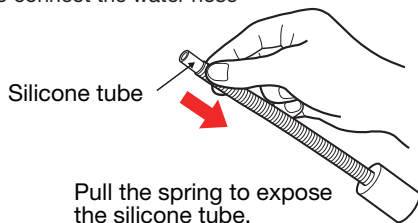Pull the spring to expose the silicone tube.

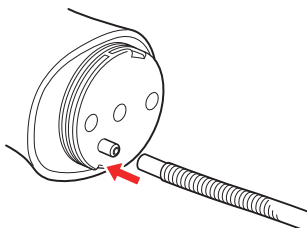
■ Storing the product

Avoid a hot, humid place such as a bathroom, and store them indoors. Otherwise, the built-in battery's life may be badly affected.

● To attach the water tank to the main unit

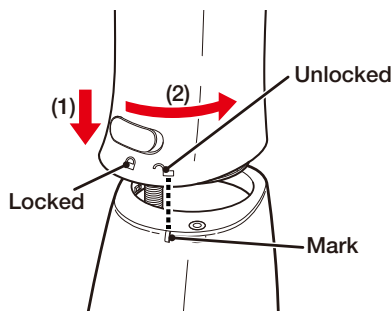To attach the water tank, place it so that the  (unlocked) mark on the main unit is aligned with the  mark on the water tank. Then, turn the main unit until the  (locked) mark on the main unit is aligned with the  mark on the water tank.

Be careful not to touch the power switch of the main unit.

Cleaning and storage (continued)

■ Lithium-ion battery

- The product has a built-in lithium-ion battery. A lithium-ion battery tends to show a decrease in performance the longer it is used, even in normal use conditions. Over time, the battery may take longer to be charged fully or may not last long as expected, but they are not a fault. The battery has reached the end of its life if it lasts only for an extremely short time when fully charged, but it cannot be replaced.

Recycling of lithium-ion batteries

Li-ion 00

- Lithium-ion batteries are a precious resource that can be recycled. We appreciate your cooperation in recycling them.
- Do not dispose of the main unit as household waste. The lithium-ion battery may crash in a garbage truck and cause a short circuit, leading to a fire or heat generation.
- Follow the guidance of the local municipal government when disposing of the main unit. For the local rules, contact your local government.
- Do not disassemble the battery when disposing of it.

Troubleshooting

● Charging

Problem	Probable cause	Remedy
Pressing the power button does not start flossing.	▶ The charging indicator is blinking red.	▶ Charge the battery. (⇒ p. 7)
	▶ The product is used outside the recommended operating temperature of 5 °C to 40 °C.	▶ Outside the recommended operating temperature, the battery protection feature may be activated, disabling the product's operation. Use it within the recommended operating temperature.
Cannot charge the battery.	▶ The plug of the USB AC adapter (commercially available) is not connected to an outlet.	▶ Connect the plug of the USB AC adapter (commercially available) to an outlet. (⇒ p. 7)
	▶ The dedicated USB charging cable is not connected to the USB port of the USB AC adapter (commercially available).	▶ Connect the dedicated USB charging cable to the USB port of the USB AC adapter (commercially available). (⇒ p. 7)
	▶ The dedicated USB charging cable is not connected to the USB charging port of the main unit.	▶ Connect the dedicated USB charging cable to the USB charging port of the main unit. (⇒ p. 7)
	▶ The dedicated USB charging cable is cut or nearly cut.	▶ Stop using the cable immediately. Purchase a USB Type-A — USB Type-C cable (commercially available) for replacement.
	▶ The product was left uncharged for a long time.	▶ A completely discharged battery can no longer be used. Charge approximately every six months to prevent it from being discharged completely.
	▶ The product is used outside the recommended operating temperature of 5 °C to 40 °C.	▶ Outside the recommended operating temperature, the battery protection feature may be activated, disabling the product's charging. Use it within the recommended operating temperature.

Troubleshooting (continued)

Problem	Probable cause	Remedy
The battery lasts only for a few minutes when fully charged.	It is considered that the battery has come to the end of its life.	The built-in battery cannot be removed or replaced. Dispose of the product by following the guidance of the local municipal government.

● During use

Problem	Probable cause	Remedy
No water jets come out, or water jets are weak.	There is not enough water in the water tank.	Pour cold water or lukewarm water (40 °C or less) before use. (⇒ p. 9)
	The flosser nozzle is clogged.	Clean the flosser nozzle. (⇒ p. 14)
	The outlet of the flosser nozzle is deformed.	Replace the flosser nozzle. (⇒ p. 9, 20)
	The water hose is deformed or damaged.	Replace the water hose. (⇒ p. 15)
	The charging indicator is blinking red.	Charge the battery. (⇒ p. 7)
Water leaks from the root of the flosser nozzle.	The flosser nozzle is wrongly attached to the main unit.	Correctly attach the flosser nozzle to the main unit. (⇒ p. 9)

*If the above issues are not resolved or if there is any other abnormality or malfunction, please contact the retailer where this product was purchased. Do not attempt to repair it by yourself.

FAQ

Question	Answer
I feel pain in soft mode. Can I continue to use it?	If you have a strong pain or intense bleeding from the gums, stop use and consult a dentist.
Can I use it in a bathroom?	Yes.
How often will I need to renew the flosser nozzle?	Replace it in approximately six months after starting use for good hygiene.
How often will I need to renew the water hose?	Replace it in approximately six months after starting use for good hygiene.
Can I use mouthwash or dental rinse in the water tank?	Do not use such a liquid to prevent malfunction.
Can I pour warm water into the water tank?	You can use lukewarm water of 40 °C or lower. Do not use water warmer than 40 °C.
Can I leave it unused for a long time?	If the battery is discharged completely, it may not work any longer. Therefore, charge it approximately every six months when leaving it unused for a long time.

Replacement parts

- Confirm the model numbers of the main unit (HC-JF03WH/HC-JF03BK) and the parts to be replaced before purchasing replacement parts.

Product name	Model number	Quantity
Dental Washer replacement flosser nozzle	HC-JF03NZ	2
Dental Washer replacement water hose	HC-JF03TB	1

[Replacement frequency]

Replace the flosser nozzle and water hose approximately every six months for good hygiene. If you find the flosser nozzle and water hose deformed (worn) or discolored within six months, replace them with new ones earlier.

Specifications

Rated voltage	5 V (USB power delivery)
Power consumption (max)	3.5 W (during charging)
Charging method	Charging via USB Type-C®
Battery type	Lithium-ion battery
Rated battery capacity	3.7 V/1600 mAh
Charging time	Approx. 3 hours (to fully charge an exhausted battery)
Number of charging cycles	Approx. 300 *The number of charging cycles is approximate. Please note that it may vary depending on the environment.
Max continuous operating time	Approx. 40 minutes (in hard mode with a fully charged battery)
External dimensions	Approx. W75 x D82 x H215 mm (excluding flosser nozzle)
Weight	Approx. 260 g (main unit only)
Tank capacity	Approx. 350 ml
Accessories	Flosser nozzle x 2 and dedicated USB charging cable (USB Type-C® — USB Type-A)
Recommended operating temperature	5 °C to 40 °C
Country of manufacture	China

★ Perform regular inspections on your aging Dental Washer.

Check for
these
symptoms
during
use.

- The main unit or dedicated USB charging cable gets abnormally hot.
- There is a strange odor or noise.
- The main unit has cracked or broken.
- You feel tingling when touching the main unit.
- Any other abnormality or malfunction

Stop
use

In case of these symptoms, disconnect the dedicated USB charging cable from the main unit to prevent malfunction or accident. Next, consult the retailer where this product was purchased.

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- ELECOM CO., LTD. owns the copyright of this manual.
- Unauthorized copying and/or reproducing of all or part of this manual is prohibited.
- If you have any concerns or question regarding this manual, please contact the retailer from whom you purchased this product.
- The images are for illustration purposes only and may differ from your actual product.
- The specifications and appearance of this product may be changed without prior notice for product improvement purposes.
- We shall bear no responsibility for any device failure or malfunction of other devices resulting from the use of this product.
- When exporting this product, check the export regulations for the country of origin.
- USB Type-C® and USB-C® are registered trademarks of USB Implementers Forum.
- Company names, product names, and other names in this manual are either trademarks or registered trademarks.

Global Support

Contact your local retailer with any product inquiries.

www.elecom.co.jp/r/s536

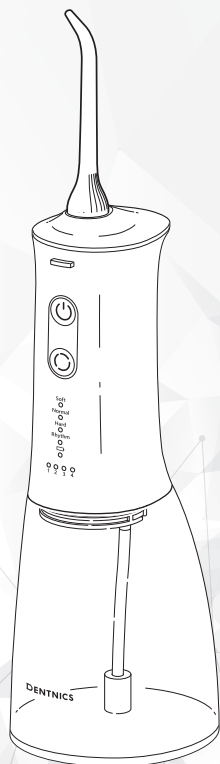
Dental Washer Tankless Compliant Model: HC-JF03
Manual, Safety precautions, and Compliance with regulations
1st edition, July, 2025
ELECOM
©2025 ELECOM Co., Ltd. All Rights Reserved.

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Trademark and name of the manufacturer: