
**Manual(English) ▶▶▶ p.2 - p.27**

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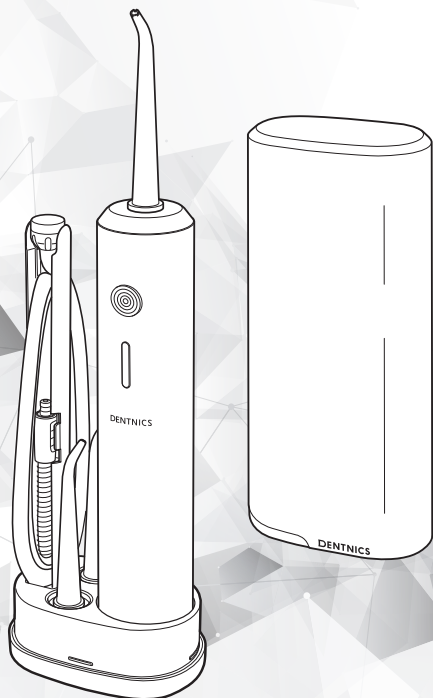
For residential use

Manual (ENGLISH)

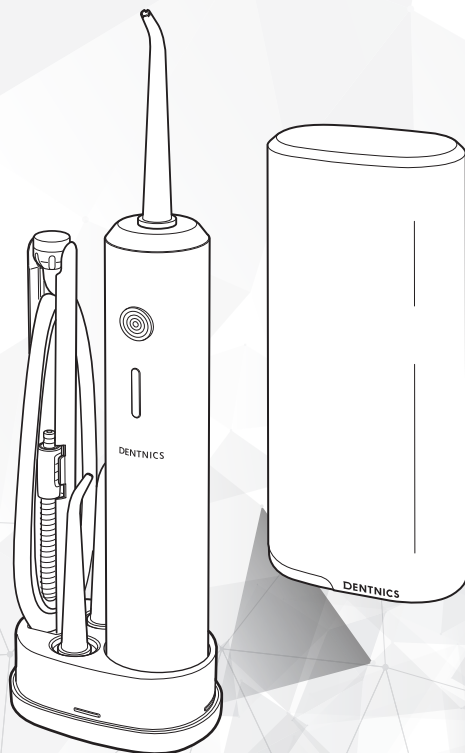
Read this manual and use the product safely and correctly.  
After reading, keep the manual in an easily accessible location for future reference.

\*In this manual, HC-JF01WH/HC-JF01BK will be referred to as “the product” except in some places.

# Package contents

The product package includes the items below. Check if they are all supplied. Although every effort has been made to ensure proper packaging, contact the retailer where this product was purchased if any shortage or breakage is discovered.

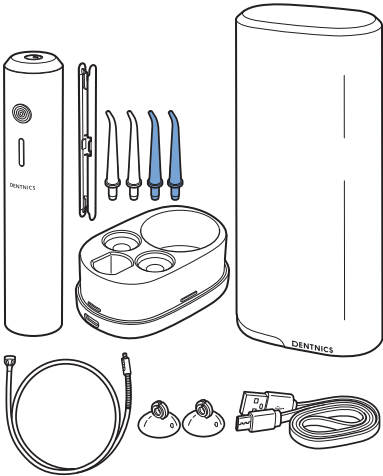- ☐ Main unit × 1
- ☐ Flosser nozzle × 4 (transparent and blue × 2 each)
- ☐ Dedicated charging stand × 1
- ☐ Cap and tank × 1
- ☐ Dedicated USB charging cable × 1
- ☐ Water hose × 1
- ☐ Water hose holder × 1
- ☐ Sucker for drying water hose × 2

## Table of contents

<b>Package contents .....</b>	<b>1</b>	<b>Cleaning and storage.....</b>	<b>14</b>
<b>Safety precautions .....</b>	<b>2</b>	After each use.....	14
<b>Names and functions of components .....</b>	<b>5</b>	Storing the product .....	15
<b>Before use .....</b>	<b>7</b>	Cleaning the product (weekly or so) .....	16
Charging the product .....	7	When the water jets are weak .....	17
<b>How to use.....</b>	<b>9</b>	Lithium-ion battery .....	18
Setting the flosser nozzle .....	9	<b>Troubleshooting.....</b>	<b>19</b>
Setting the water hose holder and		<b>FAQ .....</b>	<b>22</b>
water hose .....	9	<b>Replacement parts.....</b>	<b>23</b>
Pouring water.....	10	<b>Specifications.....</b>	<b>23</b>
Selecting a mode .....	11	<b>Regular inspections .....</b>	<b>24</b>
Flossing your mouth.....	11	<b>Compliance status .....</b>	<b>25</b>

# Safety precautions

Read the safety precautions carefully and use the product correctly. The precautions in this section describe the actions and behaviors to take to prevent hurting yourself and others or damaging property. The meaning of the symbols used in this manual is explained below.

## ● Meaning of symbols

 **WARNING** indicates matters that have a “probability of death or serious injury can be anticipated”.

 **CAUTION** indicates matters that have a “probability of injury or property damage can be anticipated”.

 General precautions

 Do not touch with wet hands

 Do not wet

 General prohibitions

 Do not disassemble

 Follow instructions

 indicates “matters that require careful attention”.

 indicates “recommendations” and “reference information”.

## WARNING

 **Do not use the product for any purpose not described in the manual.**  
Doing so may cause a fire, ignition, electric shock, burn, or injury.

 **Do not disassemble, repair, or modify the product.**  
Doing so may cause a malfunction, fire, electric shock, or injury. In the event of a malfunction, stop using the product and disconnect the dedicated USB charging cable. Next, consult the retailer where this product was purchased.

 **Do not put the product in a fire, heat it, or leave it in an extremely hot place.**  
The product contains a rechargeable battery. Being extremely hot, it may cause heat generation, ignition, rupture, or injury.

 **Do not use the product near medical or home therapy equipment.**  
Do not use the product near implanted medical devices (e.g., a cardiac pacemaker and a defibrillator), life support and management devices (e.g., a heart-lung machine and a respirator), or wearable medical devices (e.g., a wearable automatic defibrillator and an ambulatory infusion pump). Doing so may cause the devices to malfunction.

 **Do not allow children or those who are incapable of, requiring assistance in, or unfamiliar with operating the product to use it alone.**  
Doing so may cause an electric shock, burn, or injury.

 **Do not allow water to spill over the dedicated charging stand, dedicated USB charging cable, or USB AC adapter (commercially available). Do not fully wash them.**  
Doing so may cause a malfunction, fire, or electric shock.

 **Do not use or store the dedicated charging stand, dedicated USB charging cable, and USB AC adapter (commercially available) in a bathroom or any other humid place.**  
Doing so may cause a fault current or short circuit, resulting in a fire or electric shock.

## Safety precautions (continued)

### **WARNING**

- 
-  **Always use the dedicated USB charging cable.**  
Failure to do so may cause a malfunction, error, fire, ignition, electric shock, or burn.
- 
-  **Do not allow conductive foreign objects (such as metal pieces, pencil lead, or foreign particles) to get inside the USB terminal.**  
They may cause a malfunction, fire, ignition, electric shock, or burn.
- 
-  **Immediately stop using the product in the event of an abnormality or malfunction.**  
Failure to do so may cause a fire, ignition, electric shock, burn, or injury.
- There is a strange odor or noise.
  - The dedicated USB charging cable or USB AC adapter (commercially available) gets abnormally hot.
  - You feel tingling when touching the main unit or dedicated charging stand.
  - The main unit or dedicated charging stand has cracked or broken.
- If you find symptoms above, immediately stop using the product and disconnect the dedicated USB charging cable. Next, consult the retailer where this product was purchased.
- 
-  **Do not use wet hands when handling the dedicated charging stand or connecting or disconnecting the dedicated USB charging cable and USB AC adapter (commercially available).**  
Doing so may cause a fire, electric shock, or burn.
- 
-  **Always use a USB AC adapter that is USB Type-A with a rated output of 5 V/1 A or higher.**  
Failure to do so may cause a malfunction, error, fire, ignition, electric shock, or burn.
- 
-  **Do not use a dedicated USB charging cable that is damaged or has a loose USB terminal.**  
Doing so may cause a fire, ignition, electric shock, or burn.
- 
-  **Do not use the dedicated USB charging cable in such a manner that it may generate heat.**  
Doing so may cause a fire, ignition, electric shock, or burn.
- Use it while coiled up
  - Modify or damage it
  - Apply excessive force to it by pulling, bending strongly, or twisting it
  - Place a heavy object on it
  - Leave it squeezed between objects or close to a high-temperature object
- 
-  **Do not wind the dedicated USB charging cable around the main unit or dedicated charging stand.**  
Doing so may cause a broken cable, fire, ignition, electric shock, or burn.
- 
-  **Do not use the product or put accessories or parts within reach of young children.**  
Doing so may cause an electric shock, injury, or malfunction.
-

## Safety precautions (continued)


### CAUTION

-  **Do not use the product on unstable or slippery surfaces.**  
Doing so may cause an injury or malfunction.
-  **Do not drop, step on, throw, or subject this product to strong forces, shocks, or vibrations.**  
Doing so may cause a broken dedicated USB charging cable or USB terminal, resulting in a fire, ignition, electric shock, or burn.
-  **Disconnect the dedicated USB charging cable and unplug the USB AC adapter (commercially available) from the outlet except while charging the product.**  
Failure to do so may cause a malfunction, fire, or electric shock.
-  **Use the product for residential use.**  
The product is designed for residential use.
-  **Do not use the product for any purpose not described in the manual, or for young children or pets.**  
Doing so may cause a malfunction, fire, or burn.
-  **Do not use the product for any purpose other than flossing your mouth.**  
Water jets may cause an injury when applied to your eyes, nose, ears, throat, under the tongue, or any other delicate part.
-  **If you are suspected of gum disease or are concerned about the health of your teeth and mouth, consult a dentist before use.**  
Failure to do so may cause injury or damage to your teeth and gums.
-  **Do not use the product if you feel numb in your mouth.**  
Doing so may cause injury or damage to your teeth and gums.
-  **Do not allow any part other than the tip of the flosser nozzle to come into contact with your teeth and gums.**  
Doing so may cause injury or damage to your teeth and gums.
-  **Do not use hot water above 40 °C.**  
Doing so may cause a burn.
-  **Do not use any liquid other than cold or lukewarm water, such as mouthwash and dental rinse.**  
Doing so may cause a malfunction.
-  **Do not share a flosser nozzle with other people including your family members.**  
Doing so may cause an infection or inflammation through saliva and mucosa.
-  **After use, clean and dry the main unit, cap and tank, flosser nozzle, and water hose.**  
Failure to do so may allow mold and bacteria to grow, which can be harmful to your health.

To dispose of the product and package, follow the local rules in your area.

# Names and functions of components

## ● Main unit
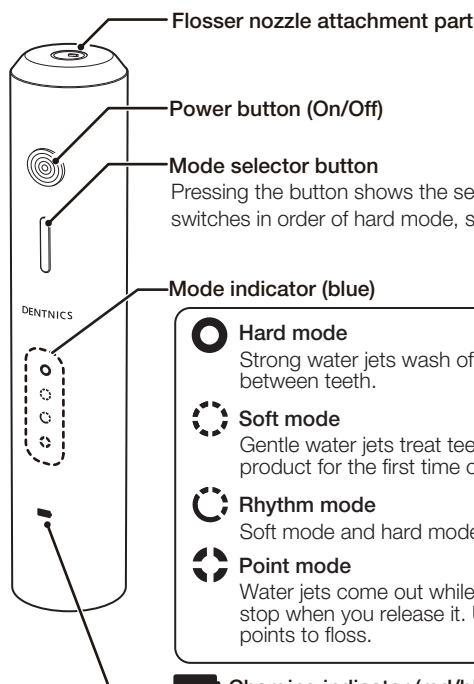
### Hard mode

Strong water jets wash off plaque on the surfaces and food particles between teeth.

### Soft mode

Gentle water jets treat teeth and gums. Optimal when you use the product for the first time or have sensitive gums.

### Rhythm mode

Soft mode and hard mode alternate to stimulate and treat gums.

### Point mode

Water jets come out while you are holding down the power button and stop when you release it. Useful when aiming water jets at specific points to floss.


### Charging indicator (red/blue)

Blinks red when the battery level is low, blinks blue when charging is in progress, and stays lit blue when charging is completed.

## ● Flosser nozzle × 4 (transparent and blue × 2 each)


You can divide them for use with your partner and family members.

\*Flosser nozzles are consumables.

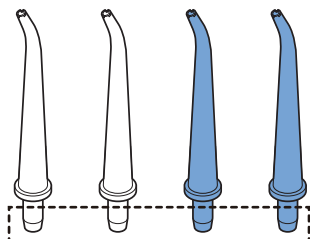Part to be inserted into the main unit

## ● Dedicated charging stand

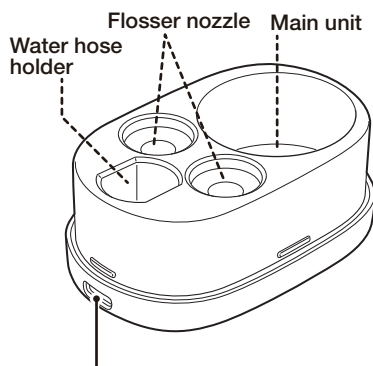
### USB charging port

Connect the dedicated USB charging cable (USB Type-C® connector) supplied.

## Names and functions of components (continued)

### ● Dedicated USB charging cable

USB Type-A connector

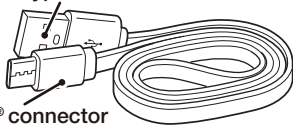
USB Type-C® connector


\*This product does not include a USB AC adapter. Use a commercially available USB AC adapter that is USB Type-A with a rated output of 5 V/1 A or higher.

### ● Water hose holder

\*A holder to store the water hose.

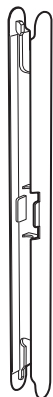
### ● Sucker for drying water hose x 2
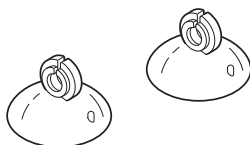
### ● Water hose

Water feed part

Connector part that connects to the main unit

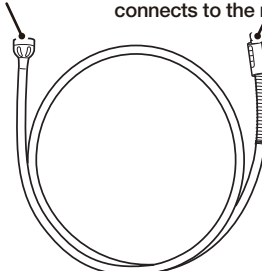
\*There may be distilled water remaining inside the main unit and water hose, which was used in the performance (water jet) test.


Rinse the product first after purchase. Refer to "[Rinsing] ⇒ p. 10" for details.

### ● Cap and tank

Use this as a tank during flossing and as a cap during storage.

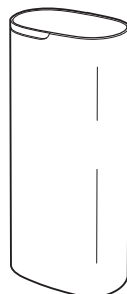When used as a tank

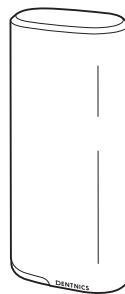When used as a cap

## Before use

### [Cautions for use]

- Do not use any liquid other than cold or lukewarm water (40 °C or lower), such as salt water, mouthwash, dental rinse, and other solutions and chemicals for mouth cleansing.
- Do not use a dryer or fan heater to dry the main unit and accessories.
- Store the product in a less humid place after use. Do not leave it in a humid place such as a bathroom.

### [Precautions before use]

#### ● About this product

- Use the product in combination with tooth brushing. Plaque building up on the tooth surfaces cannot be removed by water jets alone, so tooth brushing is essential. It is recommended that you brush your teeth before flossing your mouth with the product.  
\*Effects may vary from person to person.

#### ● Bleeding from gums

- With unusual stimulus on gums, you may experience bleeding for the first few weeks after starting use, even if you have healthy gums. If bleeding is intense or continues for longer than two weeks, it may be a symptom of abnormal gums (such as gum disease). If so, stop use immediately and consult a dentist.

#### ● Precautions about waterproofness

- The main unit is waterproof (IPX7\*) enough to allow use in a bathroom, but not completely waterproof. Do not use it under water.  
\*The product has passed the watertight test that evaluates its capability to prevent water entry that may cause adverse effects on the main unit by gradually immersing it one meter under water and leaving it for 30 minutes.
- The dedicated charging stand and dedicated USB charging cable are not waterproof. Do not charge the product in a bathroom, near a humidifier, or in any other place that may be exposed to water.

## ■ Charging the product


- If you are using the product for the first time (just after purchase) or have left it unused for more than six months and the battery becomes low during use and the charging indicator blinks red, be sure to charge it fully up to maximum.
- This product does not include a USB AC adapter. Use a commercially available USB AC adapter that is USB Type-A with a rated output of 5 V/1 A or higher.

### 1

Connect the dedicated USB charging cable to the USB charging port of the dedicated charging stand.

### 2

Connect the dedicated USB charging cable to a USB AC adapter (commercially available), and then to an outlet.

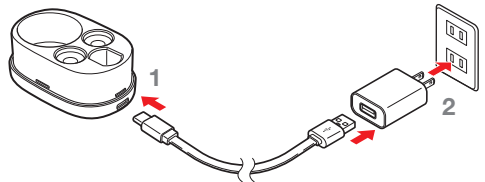
Wipe the USB charging port with a cloth and dry it well if it is wet. Leaving it wet or damp may cause a fault current.

## Before use (continued)

### 3

Set the main unit on the dedicated charging stand. The charging indicator blinks blue and charging starts.

Align the protrusion of the dedicated charging stand with the water port hole on the bottom of the main unit.
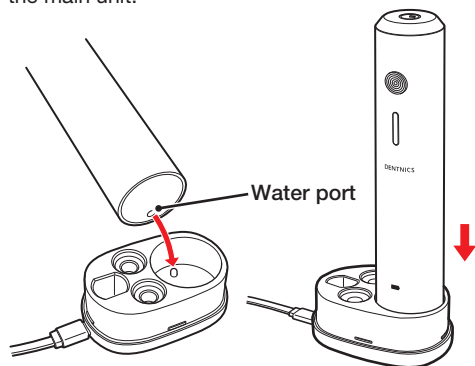
### 4

If the charging indicator stays lit blue and charging is completed, disconnect the dedicated USB charging cable from the dedicated charging stand.

- You can continue using the product for a while after the charging indicator starts blinking red, but the power will go out soon. Charging it early after blinking starts is recommended.
- If you are using the product for the first time (just after purchase) or have left it unused for more than six months, the charging time may change.
- Repeatedly charging the battery gradually shortens the operable time after a full charge. If the operable time has become extremely short, the battery has reached the end of its life. However, the battery cannot be replaced. Refer to “[Lithium-ion battery](#) ⇒ p. 18” for details.
- The product may not be charged or take more time to be charged outside the recommended operating temperature (5 °C to 40 °C).

- Immediately stop using the dedicated USB charging cable if it is broken. Purchase a USB Type-A — USB Type-C cable (commercially available) for replacement.
- Wipe water off regularly to prevent water accumulation to keep the dedicated charging stand hygienic.

#### [Charging indicator status]

Charging in progress:

The indicator blinks blue.

- It will take approximately six hours to fully charge a completely exhausted battery.

\*The product cannot be used during charging.

Charging completed:

The indicator stays lit blue.

- With the battery fully charged, the maximum continuous operating time is approximately 35 minutes (in hard mode).

\*It depends on the ambient temperature and the use conditions.

Low battery: The indicator blinks red.

- The lamp blinking red when the power is turned on means it is time to charge.

# How to use


Check the flosser nozzle and water hose for damage and deformation before use.

## ■ Setting the flosser nozzle

Hold the main unit and attach the flosser nozzle by turning it.

- Set the flosser nozzle securely so that its tip faces the front.
- The tip comes to face the front when the flosser nozzle is inserted into the water port groove with its tip facing right.

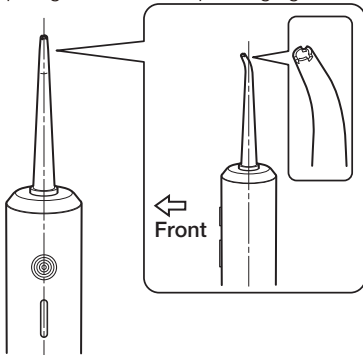
## ■ Setting the water hose holder and water hose

Align the grooves of the water feed part at the tip of the water hose with the protrusions on the inner side of the tip of the water hose holder, and gently push in the water hose to engage them.

- At the time of purchase, you will find the water hose folded on the water hose holder. Unfold and set it.
- Do not hold and pull the hose part strongly. Doing so may damage it.

Grooves of the water feed part
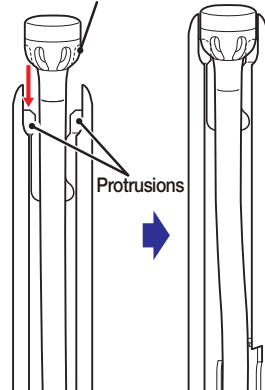
## How to use (continued)


### ■ Pouring water

1

Insert the connector that connects to the main unit on the opposite end of the water hose into the water port on the bottom of the main unit.

Make sure that it is inserted tightly without a gap.

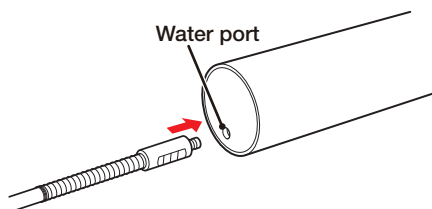
2

Pour water into the cap and tank.


- Place the cap and tank on a level surface.
- You can use a cup or plastic bottle instead of the cap and tank.

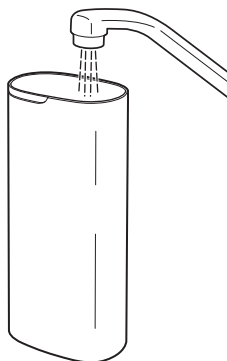- Always use fresh tap water or lukewarm water.
- If the water is too cold for you, use lukewarm water of 40 °C or lower.

3


Put the water hose holder with the water hose set to it into the cap and tank.

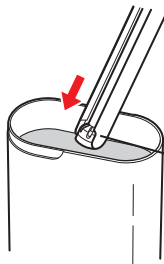Be sure to keep the tip of the water hose down and insert it.

#### [Rinsing]

- If you use the product for the first time (just after purchase) or have left it unused for more than six months, rinse it in hard mode for approximately 10 seconds.

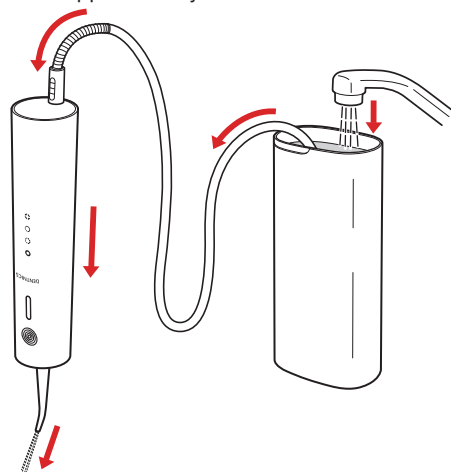The recommended operating temperature is 5 °C to 40 °C. The product may not work outside the recommended operating temperature.


## How to use (continued)

### ■ Selecting a mode

Press the mode selector button to select a mode.

- Pressing the mode selector button causes the mode indicator of the selected mode to light blue for approximately two seconds.
- To switch the mode, press the mode selector button while the mode selector button is lighting blue.
- To see the selected mode later, press the mode selector button again. The indicator of the selected mode lights blue for approximately two seconds.
- Refer to “Names and functions of components ⇒ p. 5” for the details of the modes.
- “Soft mode” is recommended when you use the product for the first time.
- The selected mode is stored after the power is turned off. The last mode used resumes the next time the power is turned on. (Memory function)
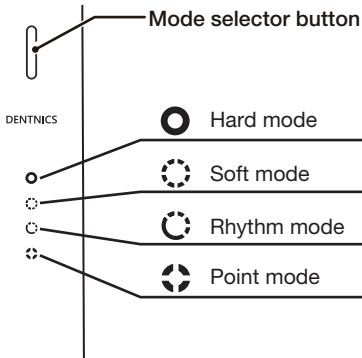
### ■ Flossing your mouth

1

Above a basin, put the flosser nozzle into your mouth.


You can also use the product in a bathroom. Refer to “[Precautions before use] ⇒ p. 7” for the waterproofness.

2

Slightly close your mouth.

Water will splash if your mouth is open. Slightly close your mouth so that water can naturally flow down into the basin.

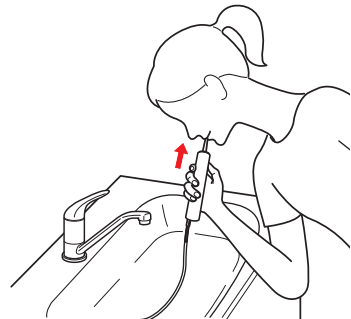
## How to use (continued)

### 3

**Press the power button to start flossing.**  
(Hold down the power button if in point mode.)

- For hard mode, soft mode, and rhythm mode, press the power button once to start and once again to stop water jets.
- If the product is not operated for approximately two minutes after the power button is pressed once or the mode is switched, water jets automatically stop.
- In point mode, water jets come out while you are holding down the power button and stop when you release it.
- You can also switch the mode by pressing the mode selector button during flossing.
- Pour water into the cap and tank or cup if it has become empty on the way.

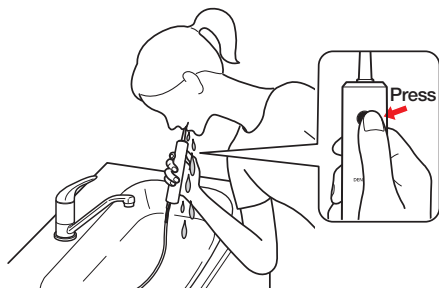
### [Tips for flossing]

- Change the mode and how to apply and move the flosser nozzle according to your purpose and the part to floss.

\*If you are concerned about your gum condition, use soft mode.


**To floss the gaps between teeth**  
(recommended mode: hard mode and point mode)

- Keep the flosser nozzle over each gap for a short time during flossing.
- Apply water jets to the upper and lower teeth from their front and back.

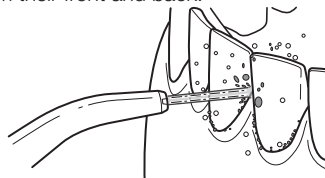
**To floss the periodontal pockets**  
(recommended mode: hard mode, soft mode, and point mode)

- Apply water jets to the gumlines.
- Slowly move the flosser nozzle along the gums.
- Floss the front and back of the upper and lower teeth and the back of the molars.

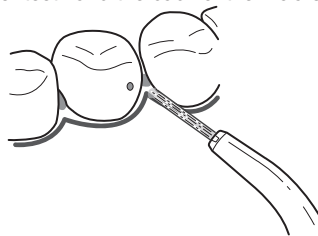**To care for the gums** (recommended mode: rhythm mode)

- Apply water jets to the gums to stimulate them.
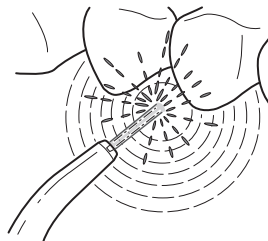
## How to use (continued)


### 4

After flossing, press the power button to turn the power off while the flosser nozzle is still in your mouth.

(Release the power button if in point mode.)

To prevent water splashing, wait until water feeding stops before taking the flosser nozzle out of your mouth.

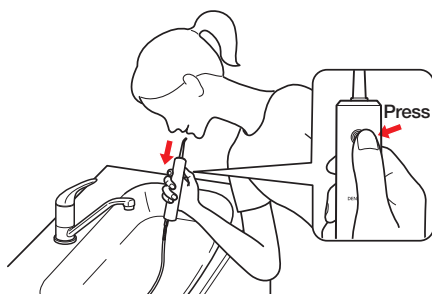If the product is not operated for approximately two minutes during continuous use, water jets automatically stop. After they automatically stop, press the power button again to resume. Keep the continuous use within approximately four minutes. After the continuous use, wait approximately 30 minutes before the next use.

# Cleaning and storage

## ■ After each use

- After each use, be sure to throw away the water and clean the product to keep it hygienic.

1

After use, rinse the product in hard mode for approximately 10 seconds.

- Refer to “[Selecting a mode](#) ⇒ p. 11” for how to switch the mode.
- Refer to “[\[Rinsing\]](#) ⇒ p. 10” for details of rinsing.

2

Remove the flosser nozzle from the main unit.


You can remove the flosser nozzle by turning it to the right and aligning it with the groove of the main unit.

3


Disconnect the water hose from the main unit.

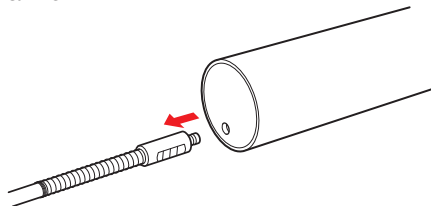Be sure to hold the connector part at the root of the water hose when disconnecting it. Pulling it by holding the hose part may damage it.

4

Remove the water feed part at the tip of the water hose from the water hose holder.
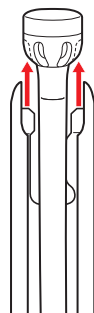
# Cleaning and storage (continued)


## ■ Storing the product

- Store the main unit and parts by standing them on the dedicated charging stand.


Avoid a hot, humid place such as a bathroom, and store them indoors. Otherwise, the built-in battery's life may be badly affected.


Store the water hose by folding it on the water hose holder.

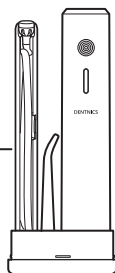- Set the water hose by pushing its tip (on the end that connects to the main unit) into the groove in the middle of the water hose holder.

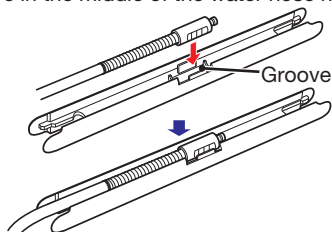
- Fold the water hose by aligning it on the upper and lower depressions of the water hose holder. Next, set the other end of the water hose by aligning the grooves.

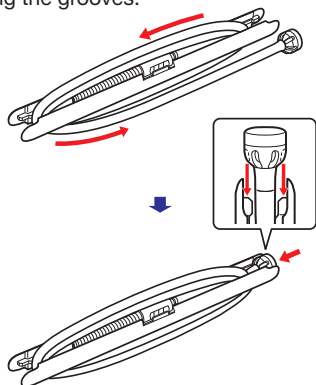
Do not pull the hose part strongly. Doing so may damage it.

- If the connector part that connects to the main unit has come off, push it in the arrow direction until it is securely in place before storing.

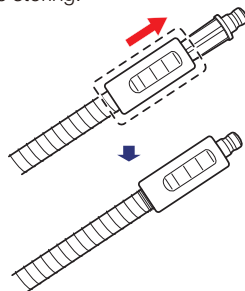- If the parts on both ends of the water hose have come off, push the hose toward the parts to reconnect them.

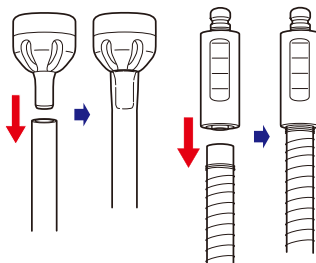- Put the cap and tank on to make it compact so that you can carry it easily.

Dry the parts completely before capping. Otherwise, water may seep out, wetting the area around it. Capping with water remaining may also result in mold and other contamination.

## Cleaning and storage (continued)

### ■ Cleaning the product (weekly or so)

- Do not use paint thinner, benzene, or alcohol. Doing so may cause a malfunction or may damage and discolor the parts.
- Do not press the power button during cleaning the product. Doing so may cause an electric shock, ignition, or malfunction.
- Do not use hot water of 50 °C or more on the main unit.
- Use neutral kitchen detergent if necessary, and rinse it off completely.


- Flosser nozzle, main unit, cap and tank, and water hose holder

1

Wash the flosser nozzle, main unit, cap and tank, and water hose holder with tap water.


If dirt is stubborn, wash them using neutral detergent diluted with water, and rinse it off completely.

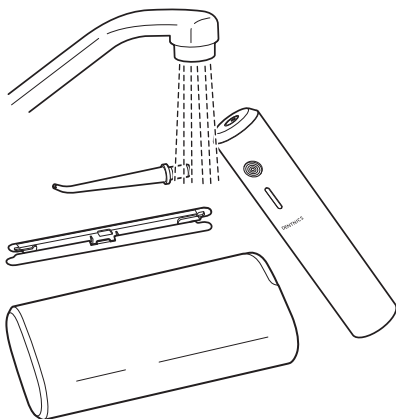Do not leave them under water.

2

Shake the main unit lightly to drain water from the inside.

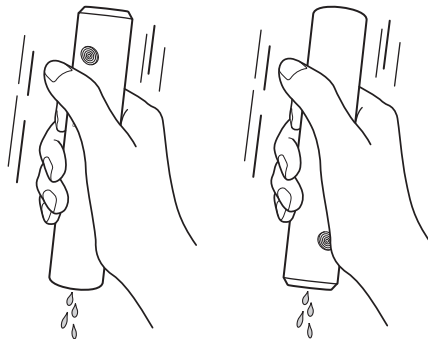
Drain from the water port on the bottom of the main unit and the flosser nozzle attachment part.

Do not hit the main unit on the basin during draining.

3

Wipe water off with a soft, dry cloth and leave the parts to dry in a well-ventilated shade.

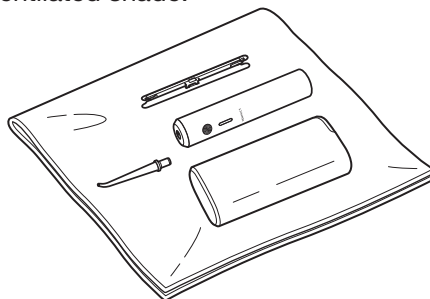
## Cleaning and storage (continued)

- Water hose

1


Wash the water hose with tap water.

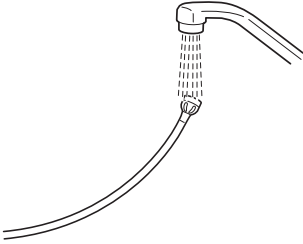2

Attach the supplied suckers for drying water hose to the wall and hang the water hose on them until it dries completely.

- Wipe dust and water off the wall before attaching the suckers.

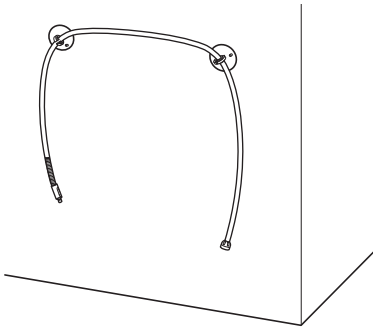- Dedicated charging stand

If it is dirty, wipe it with a soft cloth soaked in diluted neutral detergent and then with a water-wet cloth before leaving it to dry completely.

Do not wash the dedicated charging stand with water.

### ■ When the water jets are weak

If foreign materials are adhering to the water feed part at the tip of the water hose, gently rub them off with a soft brush.

## Cleaning and storage (continued)

### ■ Lithium-ion battery

- The product has a built-in lithium-ion battery. A lithium-ion battery tends to show a decrease in performance the longer it is used, even in normal use conditions. Over time, the battery may take longer to be charged fully or may not last long as expected, but they are not a fault. The battery has reached the end of its life if it lasts only for an extremely short time when fully charged, but it cannot be replaced.

### Recycling of lithium-ion batteries

**Li-ion 00**

- Lithium-ion batteries are a precious resource that can be recycled. We appreciate your cooperation in recycling them.
- Do not dispose of the main unit as household waste. The lithium-ion battery may crash in a garbage truck and cause a short circuit, leading to a fire or heat generation.
- Follow the guidance of the local municipal government when disposing of the main unit. For the local rules, contact your local government.
- Do not disassemble the battery when disposing of it.

# Troubleshooting

## • Charging

Problem	Probable cause	Remedy
Pressing the power button does not start flossing.	▶ The charging indicator is blinking red.	▶ Charge the battery. (⇒ p. 7)
	▶ The product is used outside the recommended operating temperature of 5 °C to 40 °C.	▶ Outside the recommended operating temperature, the battery protection feature may be activated, disabling the product's operation. Use it within the recommended operating temperature.
Cannot charge the battery.	▶ The plug of the USB AC adapter (commercially available) is not connected to an outlet.	▶ Connect the plug of the USB AC adapter (commercially available) to an outlet. (⇒ p. 7)
	▶ The dedicated USB charging cable is not connected to the USB port of the USB AC adapter (commercially available).	▶ Connect the dedicated USB charging cable to the USB port of the USB AC adapter (commercially available). (⇒ p. 7)
	▶ The dedicated USB charging cable is not connected to the USB port of the dedicated charging stand.	▶ Connect the dedicated USB charging cable to the USB port of the dedicated charging stand. (⇒ p. 7)
	▶ The main unit is wrongly set on the dedicated charging stand.	▶ Set the main unit on the dedicated charging stand correctly. (⇒ p. 8)
	▶ The dedicated USB charging cable is cut or nearly cut.	▶ Stop using the cable immediately. Purchase a USB Type-A — USB Type-C cable (commercially available) for replacement.
	▶ The product was left uncharged for a long time.	▶ A completely discharged battery can no longer be used. Charge approximately every six months to prevent it from being discharged completely.
	▶ The product is used outside the recommended operating temperature of 5 °C to 40 °C.	▶ Outside the recommended operating temperature, the battery protection feature may be activated, disabling the product's charging. Use it within the recommended operating temperature.

## Troubleshooting (continued)

Problem	Probable cause	Remedy
Cannot be used when fully charged.	The internal firmware may be the cause of the malfunction.	Hold down the mode selector button and the power button for six seconds. All mode indicators blink twice and the main unit is reset to the factory configuration. Next, press the power button. This may resolve the problem.
The battery lasts only for a few minutes when fully charged.	It is considered that the battery has come to the end of its life.	The built-in battery cannot be removed or replaced. Dispose of the product by following the guidance of the local municipal government.

### • During use

Problem	Probable cause	Remedy
No water jets come out, or water jets are weak.	There is not enough water in the cap and tank (or a cup).	Pour cold water or lukewarm water (40 °C or less) before use. (⇒ p. 10)
	The flosser nozzle is clogged.	Clean the flosser nozzle. (⇒ p. 16)
	The outlet of the flosser nozzle is deformed.	Replace the flosser nozzle. (⇒ p. 9, 23)
	The water hose is kinked or squeezed.	Recover the water hose to the original state.
	The water hose is deformed or damaged.	Replace the water hose. (⇒ p. 9, 10, 23)
	The charging indicator is blinking red.	Charge the battery. (⇒ p. 7)
Water leaks from the root of the flosser nozzle.	The flosser nozzle is wrongly attached to the main unit.	Correctly attach the flosser nozzle to the main unit. (⇒ p. 9)

## Troubleshooting (continued)

Problem	Probable cause	Remedy
Water leaks from the connection between the water hose and the main unit. ▶	The water hose is wrongly connected to the main unit. ▶	Correctly connect the water hose to the main unit. (⇒ p. 9, 10)

\*If the above issues are not resolved or if there is any other abnormality or malfunction, please contact the retailer where this product was purchased. Do not attempt to repair it by yourself.

## FAQ

Question	Answer
I feel pain in soft mode. Can I continue to use it?	If you have a strong pain or intense bleeding from the gums, stop use and consult a dentist.
Can I use it in a bathroom?	Yes. However, note that the dedicated charging stand is not waterproof.
How often will I need to renew the flosser nozzle?	Replace it in approximately six months after starting use for good hygiene.
How often will I need to renew the water hose?	Replace it in approximately six months after starting use for good hygiene.
Can I use mouthwash or dental rinse in the supplied cap and tank or a cup?	Do not use such a liquid to prevent malfunction.
Can I use warm water in the supplied cap and tank or a cup?	You can use lukewarm water of 40 °C or lower. Do not use water warmer than 40 °C.
Can I leave it unused for a long time?	If the battery is discharged completely, it may not work any longer. Therefore, charge it approximately every six months when leaving it unused for a long time.

## Replacement parts

- Confirm the model numbers of the main unit (HC-JF01WH/HC-JF01BK) and the parts to be replaced before purchasing replacement parts.

Product name	Model number	Quantity
Dental Washer replacement flosser nozzle	HC-JF01NZ	4 (transparent and blue × 2 each)
Dental Washer replacement water hose	HC-JF01TB	1 *A water hose holder is not included.

### [Replacement frequency]

Replace the flosser nozzle and water hose approximately every six months for good hygiene. If you find the flosser nozzle and water hose deformed (worn) or discolored within six months, replace them with new ones earlier.

## Specifications

Rated voltage	5 V (USB power delivery)
Power consumption (max)	3 W (during charging)
Charging method	Wireless charging
Battery type	Lithium-ion battery
Rated battery capacity	3.7 V/1100 mAh
Charging time	Approx. 6 hours (to fully charge an exhausted battery)
Number of charging cycles	Approx. 300 *The number of charging cycles is approximate. Please note that it may vary depending on the environment.
Max continuous operating time	Approx. 35 minutes (in hard mode with a fully charged battery)
External dimensions	Main unit: Approx. Dia.35 x H178 mm When stored (with the cap and tank set on): Approx. W76 x D48 x H190 mm
Weight	Approx. 170 g (main unit only)
Accessories	Flosser nozzle x 4 (transparent and blue x 2 each), water hose (Approx. 68 cm long), sucker for drying water hose x 2, water hose holder, dedicated charging stand, cap and tank, dedicated USB charging cable (USB Type-C® — USB Type-A)
Recommended operating temperature	5 °C to 40 °C
Country of manufacture	China

## ★ Perform regular inspections on your aging Dental Washer.

Check for  
these  
symptoms  
during  
use.

- The main unit, dedicated charging stand, or dedicated USB charging cable gets abnormally hot.
- There is a strange odor or noise.
- The main unit or dedicated charging stand has cracked or broken.
- You feel tingling when touching the main unit or dedicated charging stand.
- Any other abnormality or malfunction

**Stop  
use**

In case of these symptoms, disconnect the dedicated USB charging cable from the dedicated charging stand to prevent malfunction or accident. Next, consult the retailer where this product was purchased.

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- ELECOM CO., LTD. owns the copyright of this manual.
- Unauthorized copying and/or reproducing of all or part of this manual is prohibited.
- If you have any concerns or question regarding this manual, please contact the retailer from whom you purchased this product.
- The images are for illustration purposes only and may differ from your actual product.
- The specifications and appearance of this product may be changed without prior notice for product improvement purposes.
- We shall bear no responsibility for any device failure or malfunction of other devices resulting from the use of this product.
- When exporting this product, check the export regulations for the country of origin.
- USB Type-C® and USB-C® are registered trademarks of USB Implementers Forum.
- Company names, product names, and other names in this manual are either trademarks or registered trademarks.

Global Support

Contact your local retailer with any product inquiries.

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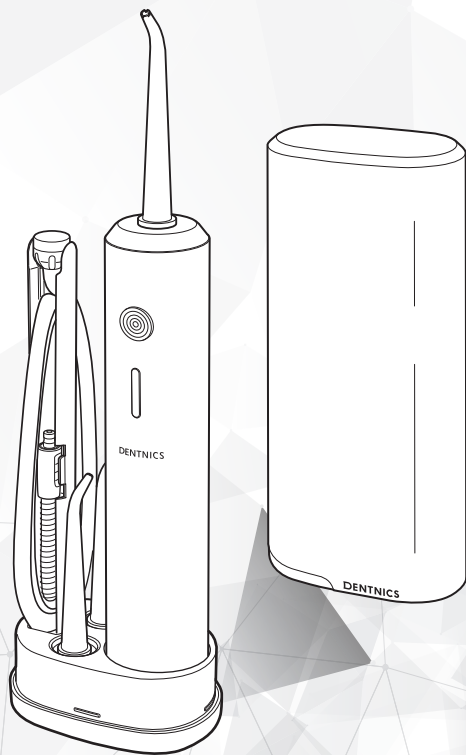
Dental Washer Tankless    Compliant Model: HC-JF01  
Manual, Safety precautions, and Compliance with regulations  
1st edition, January, 2025  
**ELECOM**  
©2025 ELECOM Co., Ltd. All Rights Reserved.

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Trademark and name of the manufacturer: